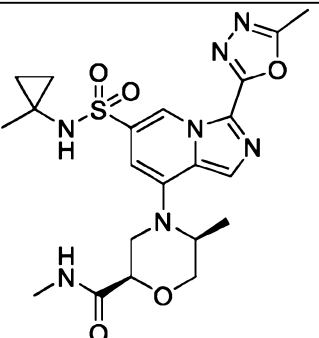
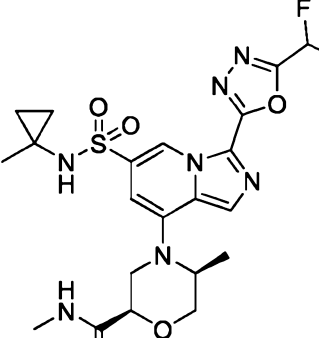
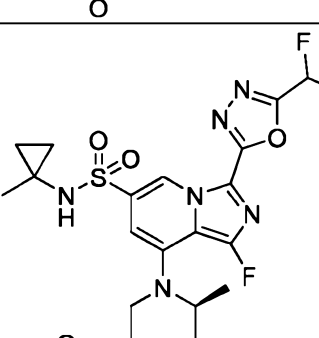
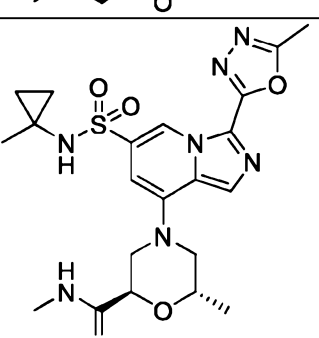


238		m/z [M+H] ⁺ 490.29
239		m/z [M+H] ⁺ 526.35
240		m/z [M+H] ⁺ 531.31
241		m/z [M+H] ⁺ 490.25

242	 <chem>CC1(C)NC(=O)C1c2cc3nc4c(ncn4C(F)F)cc3cc2S(=O)(=O)NC1CC1</chem>	m/z [M+H] ⁺ 526.35
243	 <chem>CC1(C)NC(=O)C1c2cc3nc4c(ncn4C(F)F)cc3cc2S(=O)(=O)NC1CC1</chem>	m/z [M+H] ⁺ 526.35
244	 <chem>CC1(C)NC(=O)C1c2cc3nc4c(ncn4C)cc3cc2S(=O)(=O)NC1CC1</chem>	m/z [M+H] ⁺ 477.15
245	 <chem>CC1(C)NC(=O)C1c2cc3nc4c(ncn4C)cc3cc2S(=O)(=O)NC1CC1</chem>	m/z [M+H] ⁺ 488.15

Enzyme reaction and Detection:

1. Transfer 200nL of 100x compound or DMSO to the appropriate wells of a 384 well white polystyrene microtitre plate (Corning Catalog#3574).
2. Transfer 10 uL of 2x final concentration of enzyme in assay buffer or assay buffer alone to the appropriate wells.
3. Centrifuge the plate at 1000 rpm for 30 seconds.
4. Incubate the plate at room temperature for 1 hour.
5. Transfer 10 uL of 2x substrate in assay buffer to all test wells.
6. Incubate the plate at room temperature for 10 minutes
7. Transfer 10 uL of 3x mixture of 42nM detection antibody and 2.25nM streptavidin-europium in 50mM Tris-HCL pH 7.4 to all test wells.
8. Incubate the plate at room temperature for 1 hour.
9. Read the plate on a plate reader (Envision)

Excitation: 317nM

Emission: 620nM

Emission: 665nM

Data Analysis:

[0641] EC₅₀ values were calculated in Collaborative Drug Discovery vault (CDD). Curves were fitted by CDD as response (%) vs compound concentration (uM) using a 4-parameter inhibition model using Formula 1.

[0642] Formula 1:

$$\text{Fit} = (A + ((B - A) / (1 + ((C/x)^D))))$$

$$\text{Res} = (y - \text{fit})$$

[0643] The TR-FRET EC₅₀ values for a representative number of compounds of Formula (I) or Formula (I') are provided in Table 2 below.

TR-FRET

**** <= 0.1 μM, 0.1 μM < *** <= 0.3 μM, 0.3 μM < ** <= 1 μM, 1 μM < * < 10 μM

Table 2

Example	TR-FRET EC ₅₀ (μM)	Example	TR-FRET EC ₅₀ (μM)	Example	TR-FRET EC ₅₀ (μM)	Example	TR-FRET EC ₅₀ (μM)
1	****	20	****	39	***	58	****
2	****	21	****	40	****	59	***
3	****	22	***	41	****	60	****
4	****	23	****	42	****	61	***
5	***	24	****	43	***	62	***
6	****	25	****	44	****	63	****
7	****	26	**	45	****	64	****
8	>10 μM	27	****	46	****	65	****
9	>10 μM	28	**	47	****	66	****
10	>10 μM	29	***	48	****	67	***
11	****	30	****	49	****	68	**
12	****	31	****	50	****	69	****
13	****	32	***	51	***	70	****
14	*	33	****	52	**	71	****
15	***	34	****	53	****	72	****
16	*	35	****	54	****	73	***
17	*	36	****	55	****	74	****
18	****	37	***	56	***	75	****
19	****	38	****	57	****	76	****
77	****	105	****	133	****	161	***
78	****	106	***	134	**	162	****
79	****	107	****	135	****	163	****
80	****	108	***	136	****	164	****
81	****	109	****	137	****	165	****
82	****	110	****	138	****	166	****
83	****	111	****	139	****	167	****
84	****	112	****	140	****	168	***
85	****	113	****	141	**	169	****